

# Experience Freedom

with The Golden Companion Scooter.

Whether you and your Companion are staying at home, visiting friends or grocery shopping, the Golden Companion gives you the freedom and mobility you demand and deserve.

Something truly special has happened at Golden Technologies.

We have added new features that improve performance, comfort and safety. We have unveiled new colors and updated styling.

We have redefined what you have asked for in a personal mobility scooter.

**Ask about the Golden Payment Plan**

# Experience a NEW TRADITION in a Value Added THE GOLDEN COMPANION SERIES

## New Improved Control Panel

Attractive, with new battery indicator, user friendly light/horn buttons, diagnostic diode and new key switch.

## Self Diagnostic Indicators

State of the art trouble shooting system for the consumer or technician.

## Sleep Mode (Auto Shut-Down)

The DS70 sleep mode de-energizes the power relay when the scooter panel is not used for 60 minutes.

## Stylish Recessed Headlight

Provides for additional safety. Illuminates more than other major manufacturers.


## Shorter Contoured Grips

More comfortable. Helps improve maneuverability.

## Increased Leg Room

Offers more leg room than any other scooter on the market. "Small to Tall - We Fit Them All!"

## Front Basket Standard

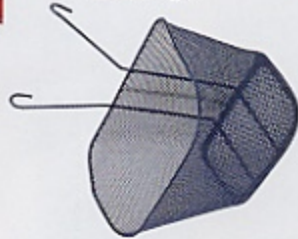

## Foam Filled Front Tire

Maintenance free. (Companion only.)

## Effortless Infinite Adjustment Tiller

New streamlined patent pending tiller allows the user to adjust/position to any comfort level without difficulty.

## Removable Automobile Style Carpet

Easy to maintain and keep clean.

## Unique Independent Front Suspension (4 Wheel Only)

Allows for a superior ride over rough and uneven terrains.

## 3 Amp Battery Charger

Enables the battery to receive a stronger/quicker charge.



## Microprocessor based, Service Friendly Electronics

Dynamic Controls DS70 motor control board offering load compensation, board diagnostics, motor and battery protection and sleep mode. *Load Compensation:* Allows the user to maintain a steady speed when traveling up hills/over varying terrain's such as gravel to side walk and carpet to dirt. *The DS70 diagnostic program:* uses flashing L.E.D. (light emitting diode) to pinpoint electrical problems if they occur. One L.E.D. is located on the of the motor board the other on the control panel.

Liftchair.com 1-888-850-5438



## New Free Wheel Brake Release

The free wheel lever, located in the rear of the unit, releases the brake and the electronics. No longer need to release to a toggle switch on the electronics tray.

**Also!**

**Increased Weight Companion 30 Companion II 3**  
Extremely wide cross-



Power Seat



Cane Holder



Rear Basket



Cup Holder



Oxygen Holder



Docking Device



Quad Cane Holder



Walker Holder



Rear View Mirror



Deluxe Seat

## Scooter Accessories and Options